

# AniMat

Double-sided veterinary heating mat

## For Safe and Effective Operation

### ► Easy to use:

**AniMat** is preset to provide the appropriate heat depending on the animal:

**AniMat** uses the animal's body temperature regulation, automatically preventing the temperature from exceeding 104°F (40°C).

### ► Convenient:

**AniMat** can be easily moved from one cage to another or from a surgery table to a cage.

The table heating mat unit allows surgeries to be performed while warming the animal in a preventive way.

### ► Secured:

**AniMat** adapts itself perfectly to the animal, due to its two different sides

- **Grey side:** for more than 20 Kg

- **Blue side:** for less than 20 Kg

Reducing risks of Ischemia and consecutive burns. The power cable is strengthened on the section accessible to the animal.

### ► Robust:

**AniMat** is made from high quality resistant materials, meeting veterinarians' needs.

The power cable is strengthened while the cover material, fully waterproof, withstands bites and scratches.

### ► Hygienic:

**AniMat** can be easily cleaned with water and soap and after drying disinfected with alcohol.

AniMat helps to maintain the animal's temperature. AniMat has been designed to be particularly robust while also providing an optimum heat distribution evenly on its entire surface.



# Veterinary Heating Mat

## ► Fields of application:

### Preventive:

Maintaining the temperature of the body is very important during and after anaesthesia. 100 % of patients get cold during anaesthesia.

Consequences of hypothermia:


- Difficulty waking up
- Risk of shock
- In the long term, risk of heart problems, renal insufficiency, etc.
- Risk of death

Awareness of hypothermia is also important during simple acts such as dental scaling. The smaller the animal is, the more affected it will be by hypothermia. It is also very important with rodents or birds.

### Curative:

Hypothermia in conjunction with trauma or anesthesia must be treated systematically. It is important to maintain the body temperature between 38 and 39°C (100 and 102°F). The heating mat enables gentle warming of the patient without any risk of burning.

Using the **AniMat** during anaesthesia or resuscitation with regular monitoring of the animal's temperature is simple but essential.

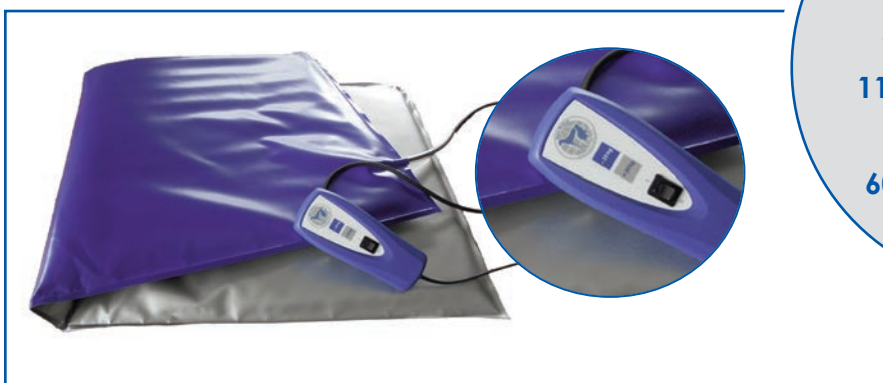

AniMat is available  
in 2 sizes:

- **Surgical table size:**  
110 x 50 cm (44" x 20")
- **Cage size:**  
60 x 50 cm (24" x 20")

### Technical features:

Surgical table size:  
Dimensions: 44" x 20" (110 x 50 cm)  
Operation: 50 W - 24 V

Cage size:  
Dimensions: 24" x 20" (60 x 50 cm)  
Operation: 27 W - 24 V



15, rue du Noroit  
22100 Taden - France

Tél. +33 (0)2 96 85 86 76

Fax +33 (0)2 96 85 86 77

[welcome@manomedical.com](mailto:welcome@manomedical.com)

[www.manomedical.com](http://www.manomedical.com)